

## How to install Tele-pointer client and setup

### 1. Download application and install

Download client software from the following URL.

<http://www.temdec.med.kyushu-u.ac.jp/data/tele-pointer-10.lzh>

Start “Setup.Exe” to install, and then follow as instructed.

Refer to see the manual which is provided by system provider.

### 2. Describe about software image

**Set your PC screen's resolution to 1280by768 before connection to adjust each site's cursor position.**

Click the icon to start software (fig. 1) then you can find dark square area at the top left corner of your desktop (fig. 2).

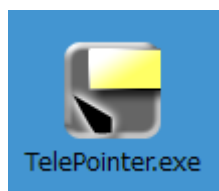


Fig. 1 Client Software

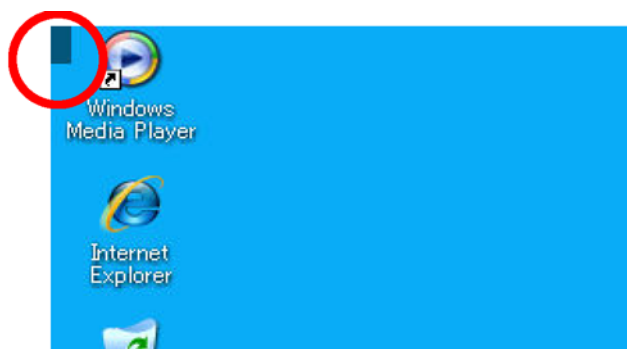






Fig. 2 Dark square area at the top left corner




### 3. Setup window shows on your desktop when you put on your cursor on the area (fig. 3).



Fig. 3 Setup window

Describe about each button

Button	Function
	Login the annotation server. This button changes “LOGOUT” when connected.
	Logout the annotation server. This button backs to “LOGOUT” after disconnected.
	Start / Finish cursor sharing. This button becomes valid during connection.
	Select appear / disappear User ID of each cursor. This button becomes valid during connection.

	Select cursor frame rate from tree options (High, normal, low). This button becomes valid during connection.
	Select setup window's display location from two options (left, right).
	Quit Client software.

#### 4. Login

Login windows appears as fig.4 when you push login button.



Fig. 4 Login Window

Input ID, Password and Server IP address correctly as informed by TEMDEC and push connect button to login.

#### <Port number settings>\*

Port number settings are changeable in this window but you do not need change anything in this window if you get some information from us.

